



FURIX

FMJQ

**Vacuum Steam Autoclave Medical
Dental Autoclave Sterilizer+Printer**

Product Features

- 01 Automatic Steam Sterilizer is a precision instrument specializing in the operation occasion in need of simple, fast and efficient sterilization by the doctor or his/her assistants, and the doctor shall take the responsibility for such work.
- 02 The sterilizer adopt microprocessor to enable intelligent control, which make easy to operate. Fully automated control, dynamic parameters and status display. Automatic fault diagnosis, over-temperature, over-pressure auto-protection to ensure the reliability of sterilization and disinfection.
- 03 The accessories of Sterilization include: Sterilization Chamber, Door unit, Insulation unit, Pipe system, Safety valve, Control system, Storage tanks, waste water tank, Steam generator and so on.
- 04 3 time pulsating vacuum make less residual air in chamber to Ensure better sterilization effect.
- 05 LCD screen show all sterilization data, which ensure better control of the sterilization process.
- 06 Safety valve installed in sterilizer, when pressure beyond 0.27Mpa. The balance pressure can be released by safety valve to ensure the operate safety.
- 06 The instrument have self-protect feature. When water tank or drainage water tank is full, the sterilization will complete pending sterilization in case the implementation of sterilization procedures, But the machine can not start to next operation unless add water or drainage. There will be buzzer to remind users to add water or drainage.
- 07 This product can be applied in the tooth device, operation device, Sanitary, Material and other heat-resisting articles existing in the Dept. of Stomatology, Ophthalmology Surgery and lab.

Parameter

- ▶ Disinfection methods: three times pulsation vacuum high temperature sterilization
- ▶ Inner chamber of the Sterilizer: 23L: 250mm X 450mm
- ▶ Voltage: 110V/60Hz, 220v 50hz both available, we will ship according to different location
- ▶ Power: 2300w
- ▶ Sterilization (Pressure/Temperature):
 - 90~130KPa/121°C (standard atmospheric pressure)
 - 20~230KPa/134°C (standard atmospheric pressure)
- ▶ Operating panel: push-button operation + LCD panel
- ▶ Panel material: high temperature matte anti-skid material
- ▶ Body material: 304 stainless steel inner liner



Volume	12L
Chamber Size	200mm X 360mm
Display	LCD
Product Size	560*450*410 mm
N.W/G.W	45/49 KG
Printer	Optional

Volume	18L
Chamber Size	250mm X 350mm
Display	LCD
Product Size	560*450*410 mm
N.W/G.W	48/50 KG
Printer	Optional

Volume	23L
Chamber Size	250mm X 450mm
Display	LCD
Product Size	560*450*410 mm
N.W/G.W	50/55 KG
Printer	Optional


PROG	Unwrapped	Unwrapped	wrapped	wrapped	Cotton	Plastic	B&D test	Vacuum test
Temperature(°C)	134	121	134	121	121	121	134	105
Pressure(KPA)	210	110	210	110	110	110	210	20
Time (min)	4-20	15-30	4-20	15-30	15-30	15-30	3.5	5

Furix Trade Venture Limited

Address: 6/F MEGA CUBE, WANG KWONG ROAD, KOWLOON, HONG KONG,
Mob: +852 9579 0472 E-mail: vet.furix@gmail.com, sales@myfurix.com
Tel: +852 9579 0472 www.myfurix.com

FURIX

Product Model	Product Volume	Product Size	Weight
WR-M/12L	200*360mm	56*45*41cm	46KG
WR-M/16L	230*350mm	56*45*41cm	48KG
WR-M/18L	250*350mm	56*45*41cm	50KG
WR-M/23L	250*450mm	68*45*41cm	55KG

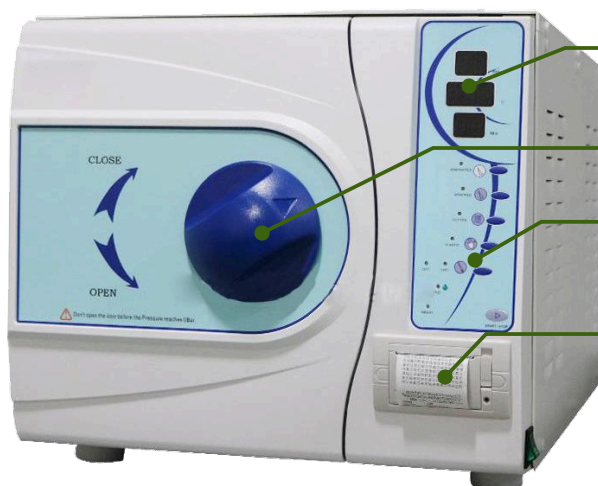



DIGITAL DISPLAY

MECHANICAL SWITCH KNOB

OPERATION BUTTON

REAL-TIME DATA PRINTER

Sterilizer built-in basket and plate size

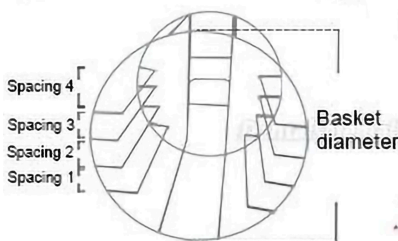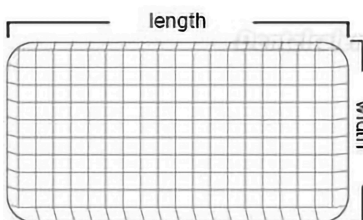

CY Size	12L	16L	18L	23L
Large plate	30.5CM*17.3CM	29.5CM*21.0CM	29.5CM*22.7CM	39.5CM*22.7CM
Medium plate	30.5CM*14.7CM	29.5CM*18.3CM	29.5CM*21.0CM	39.5CM*21.0CM
Small plate	30.5CM*14.7CM	29.5CM*18.3CM	29.5CM*20.0CM	39.5CM*20.0CM
Basket diameter	18.5CM	21.8CM	23.8CM	23.8CM
Spacing 1	3.2CM	4.4CM	2.9CM	2.9CM
Spacing 2	4.3CM	4.7CM	2.9CM	2.9CM
Spacing 3			2.9CM	2.9CM
Spacing 4	3.1CM	4CM	5CM	5.3CM

The above dimensions are measured by hand. Due to different personal measurement methods, there may be some errors, which are for reference only.

Furix Trade Venture Limited

Address: 6/F MEGA CUBE, WANG KWONG ROAD, KOWLOON, HONG KONG,
 Mob: +852 9579 0472 E-mail: vet.furix@gmail.com, sales@myfurix.com
 Tel: +852 9579 0472 www.myfurix.com

FURIX