



# FURIX

Model : FVA-5

## Furix Animal Anesthesia Machine

## 1 Top Lid

Position of Monitor (10kg MAX)

## 2 O<sub>2</sub> Flush

High Oxygen Flow (25-75L/min)

## 3 Accurate Flowmeter

Built-in Pressure Stabilizer for steady Oxygen Supply

## 4 Anesthesia Vaporizer

Accurate Pressure, Flow and Temperature Compensation  
High Capacity (300mL)  
Maintenance Free



## 5 Highly Integrated Breathing Circuit

Prevent Air Leakage and Reduce Waste of Oxygen and Anesthetic Drugs.

### a APL Valve(2-70cmH<sub>2</sub>O)

Precise Airway Pressure Control Easy to Adjust

### b ACGO Switch-off Valve

Open and Close Breathing Circuit with a Switch  
Reduce Risk of Leakage and Loop Pollution

### c Breathing Bag Connector

0-90 Degree Adjustable Rebreathing Bag Connector, which Cater to Doctor's Habits

### d Leakage Test

For Leakage Test and to Prevent Cross-infection Caused by Hand Blocking Testing

### e CO<sub>2</sub> Absorbing


Ergonomic Design with Handle for Single Hand Disassembly  
Easy to Operate

More Thorough in CO<sub>2</sub> Absorption and more Economic

### f Sensitive Respiratory Valve Diaphragm

Ceramic Material Breathing Gas Valve Sheet, Sensitive in Feedback,  
More Suitable for Low Flow Anesthesia

## Accessories

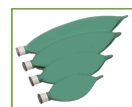

Breathing Bag  
(0.5L, 1L, 2L, 3L)



Respiratory Circuit Tube



Basket



Veterinary Anesthesia  
Mask



Oxygen Hose

## Main Specifications

- Applicability: 1kg to 100kg Animal
- Operation Mode: Semi-closed, Open Breathing Circuit for Option
- Machine with Trolley: 60\*60\*130(cm) Main Machine: 43\*23\*45(cm)
- Plimit(APL): 0.2-7.0Kpa , Flush: 25L/min-75L/min , Flow Range: 0-4L/min
- Ventilation Modes: Autonomous Breathing, Manual
- Tidal Volume: Autonomous Breathing: According to the Strength of Animal Intuitive Respiration
- Controlled Manually: 10ml-2000ml
- CO<sub>2</sub> Absorber: 1.2L, Circulatory Anesthesia is Available
- Vaporizer: Sensitive Flow, Pressure and Temperature Compensation High Quality Materials for Superb Performance, Especially at Low Flow Isoflurane, Sevoflurane, Enflurane, Halothane Tec 7 Selectatec & Tec 3 Cagemount Mount



**Shenzhen Falcon and Furix Lifesciences Co. Ltd**

**Furix Trade Venture Limited (HK)**

Address: 6/F MEGA CUBE, WANG KWONG ROAD, KOWLOON, HONG KONG,

Mob: +852 9579 0472 E-mail: sales@myfurix.com

www.myfurix.com

# FURIX