




FURIX

FVP-30

Furix Veterinary Infusion Pump

Product Specification

Flow Rate range	0.1-2000.00ml/h, increment 0.01ml/h.	Night mode	setting and adjustable.
VTBI	0.1-9999.99ml, Time: 1s-99h59min.	Temperature	-20°C-55°C
KVO	0.1-5.0ml/h, increment 0.1ml/h	Relative humidity	10%~95%
Accuracy	±5% (use the IV set calibrated)	Atmospheric pressure	40-120Kpa
Purge	0.1-2000ml/h	IV set	Dia3-5mm or setting
Infusion Mode	Rate Mode/Time Mode/Weight Mode/Sequential Mode	Bolus Rate	0.1-2000ml/h, Automatic/Manual bolus available
IV set	2 built-in brand for choice	Anti-bolus	Automatically reduce the tube pressure when the tube is blocked.
4 user defined	Common IV set / Specified IV set/ Light-proof IV set/ Precision IV set	Bubble Size	Adjustable 50/100/250/500/800 ul. Ultrasonic detection
History Record	2000 above	Occlusion Pressure	Adjustable Low/Medium /High.Double Pressure detection
Alarm	Air in Line / Downstream Occlusion/UPstream Occlusion / No Battery/VTBI Completed/KVO Finished/System Error / No Infusion Tube / Standby Time Expired / Reminder / Low Battery /No Battery Inserted/Time Near End / No AC Power.		

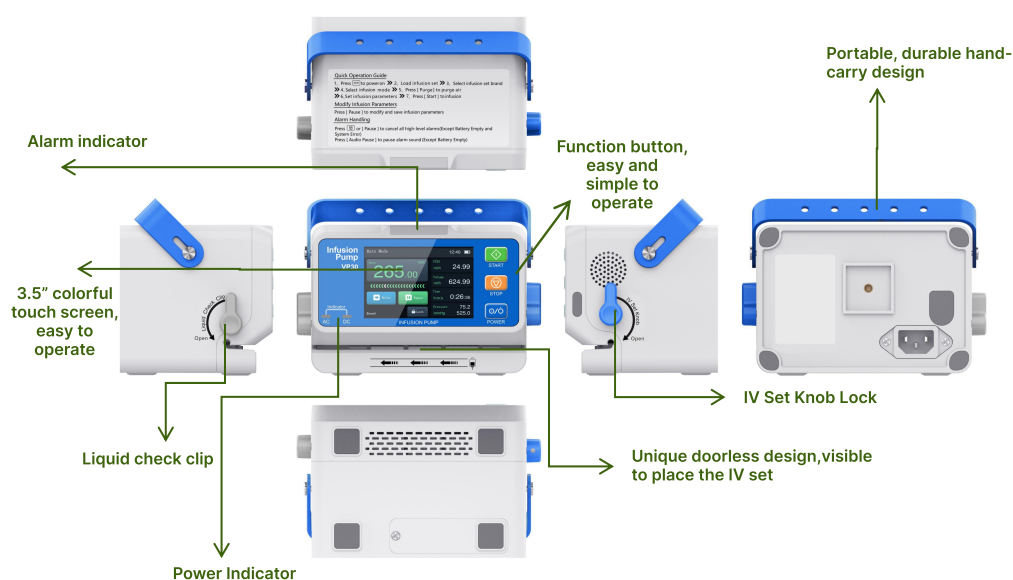

Locker design suitable for hanging on the infusion IV pole or pet cage.

- ◎ 3.5" colorful touch screen, easy to operate with human-computer interface.
- ◎ Unique doorless design, visible place the IV set directly, convenient and quick. Dual CPU control.
- ◎ IP34 dustproof and waterproof design. AC/DC power, Larger battery, long working time.
- ◎ Reduce the bubbles error alarm because horizontal place the IV set design.
- ◎ Very Light weight only 1.3 KG, small and exquisite, convenient to carry and transport.
- ◎ DPS dynamic pressure monitoring, real-time display of the current pressure value.
- ◎ Assembled German high-end medical grade motor driver chip. The unique patented silent technology reduces noise to 12dB, keeping stable operation at high and low speed, anti-shake, Power Saving, Micro-stepping and more accurate infusion, Reduces resonance, Saves energy and reduces consumption, and Efficiently cools down to keep the motor at a constant temperature.
- ◎ New independent 3 stage peristaltic device design with exclusive patented technology and First one in the global market, which it making infusion set stress evenly and improving infusion accuracy and extending the life of infusion pump and IV set, reduce workload and make the patient more comfortable in infusion process.

Furix Trade Venture Limited

Address: 6/F MEGA CUBE, WANG KWONG ROAD, KOWLOON, HONG KONG,
 Mob: +852 9579 0472 E-mail: vet.furix@gmail.com, sales@myfurix.com
 Tel: +852 9579 0472

FURIX